

# Separating two signal based on dcr (detector channel ratio) in paraview


## Goal:

- to separate signal from 2 detectors (CY5 near vs CY5 far), especially when you have two dyes for example AF640 and AF680



1. Open image in **Inspector**

2. **Minflux Data Panel** → Data source: it is usually **None**, so you must select data source that you want to view in Paraview.


3. **Minflux Data Panel** → **Launch Paraview**

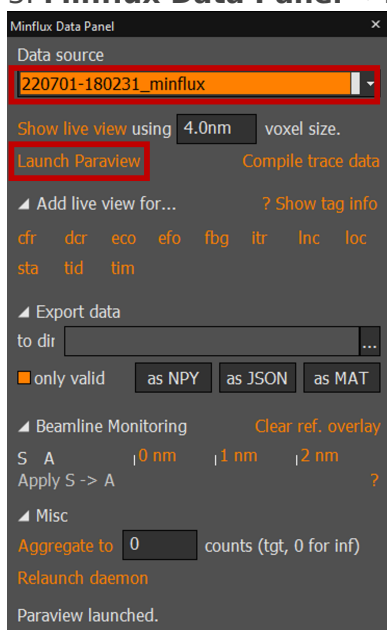



4. in **Paraview** → under localization tab → select coloring to dcr (detector channel ratio)

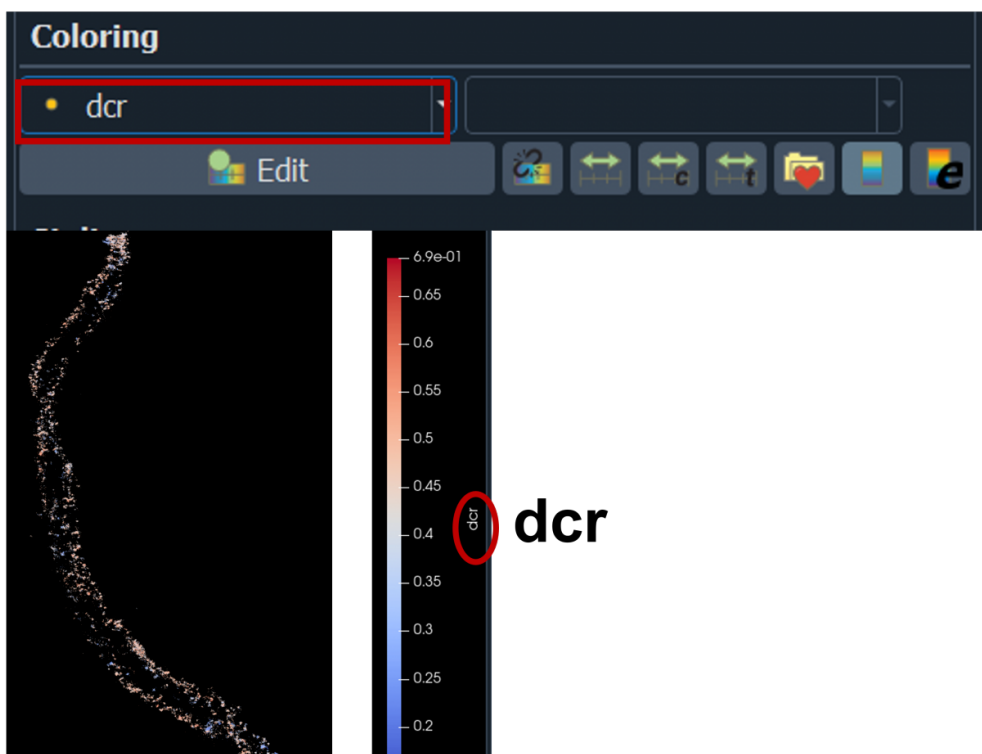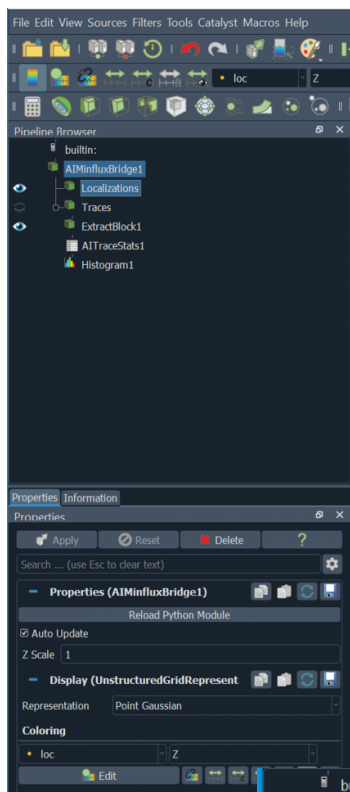


5. in **Paraview** →  
input array dcr →

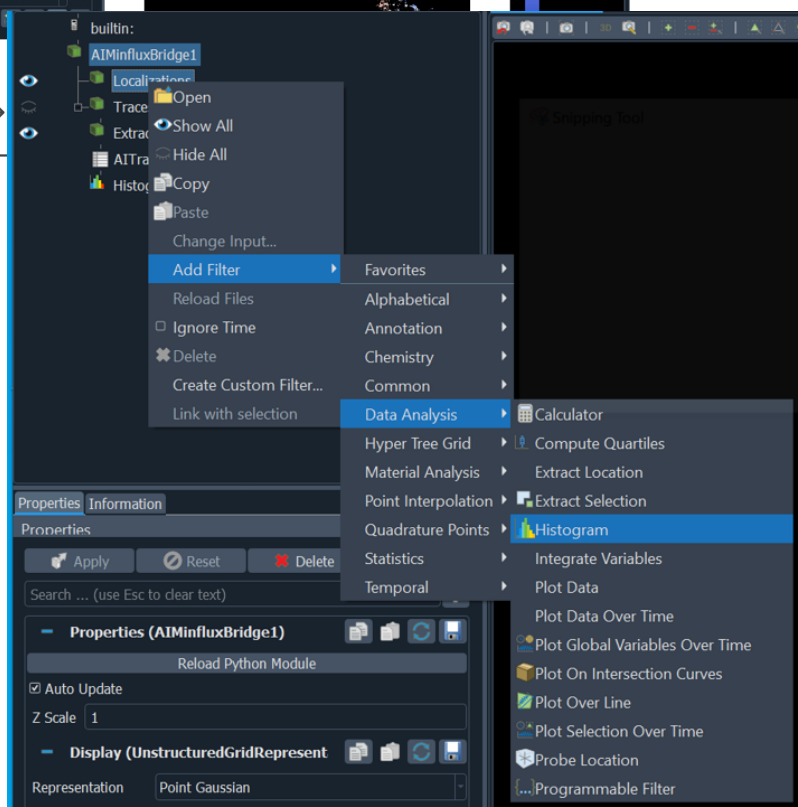


analysis → histogram → **select**

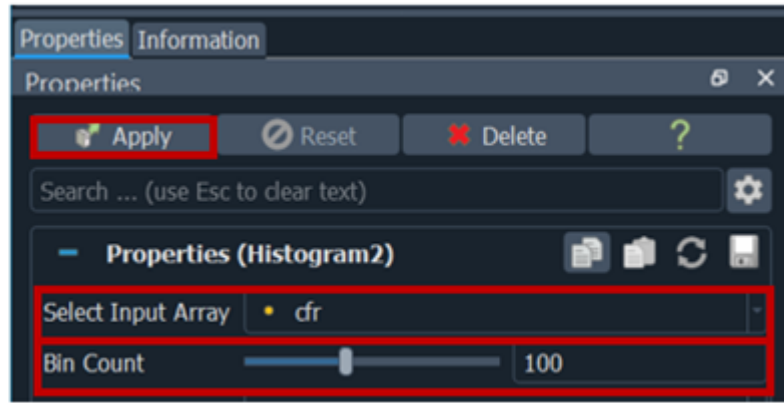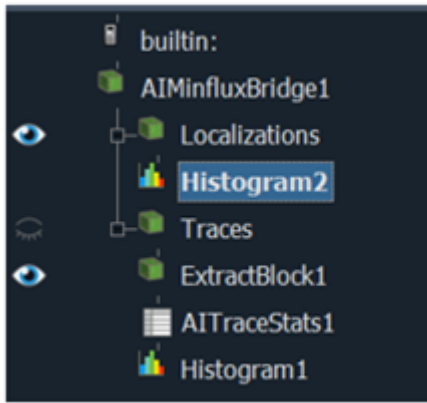

6. From histogram plot we can see two peaks from CY near and CY far detector, we can separate the two channels for example here we want to discard dcr value from 0.32 to 0.36

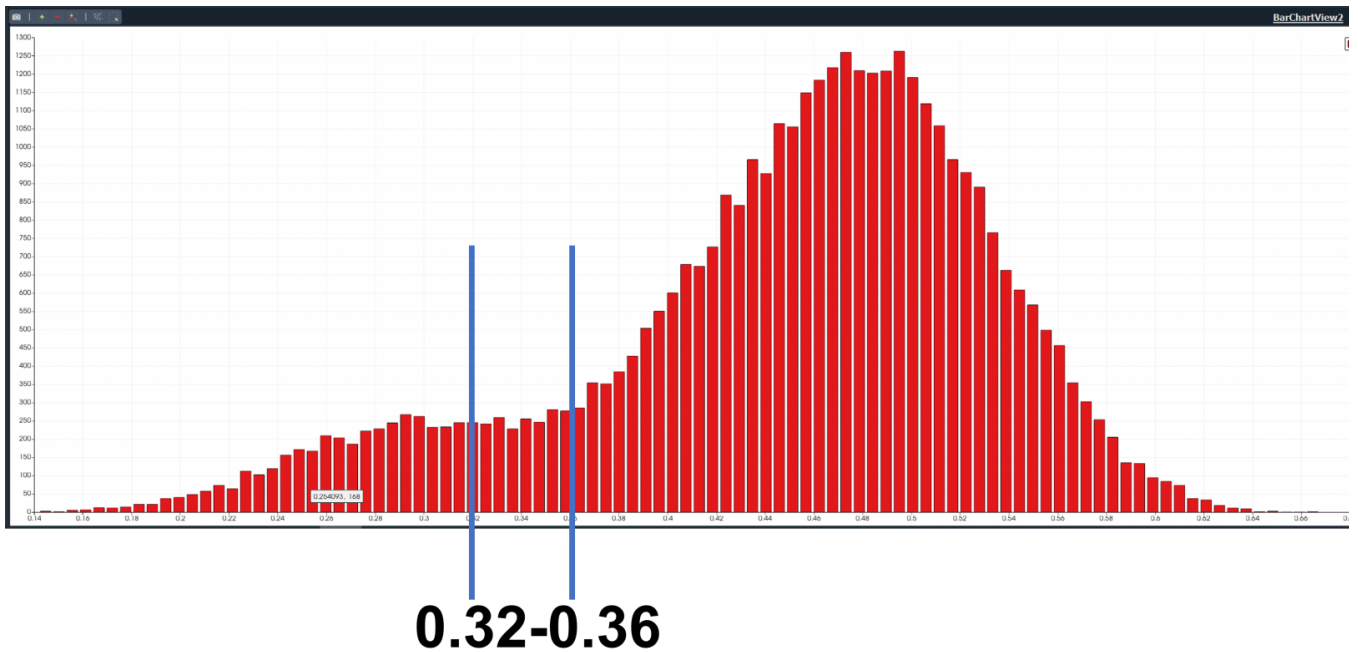


7. localizations tab: you will see color map editor or rainbow icon (edit color map) to show the color map editor

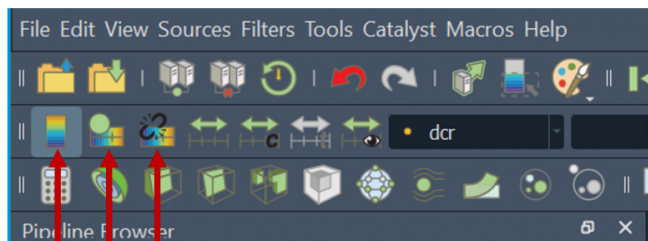


Use Separate Color Map  
Edit Color Map


Toggle Color Legend

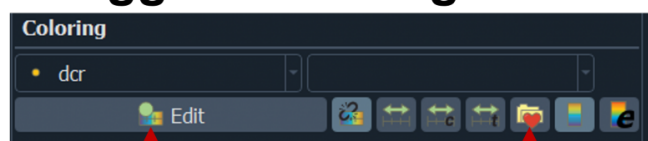

Edit Color Map Choose

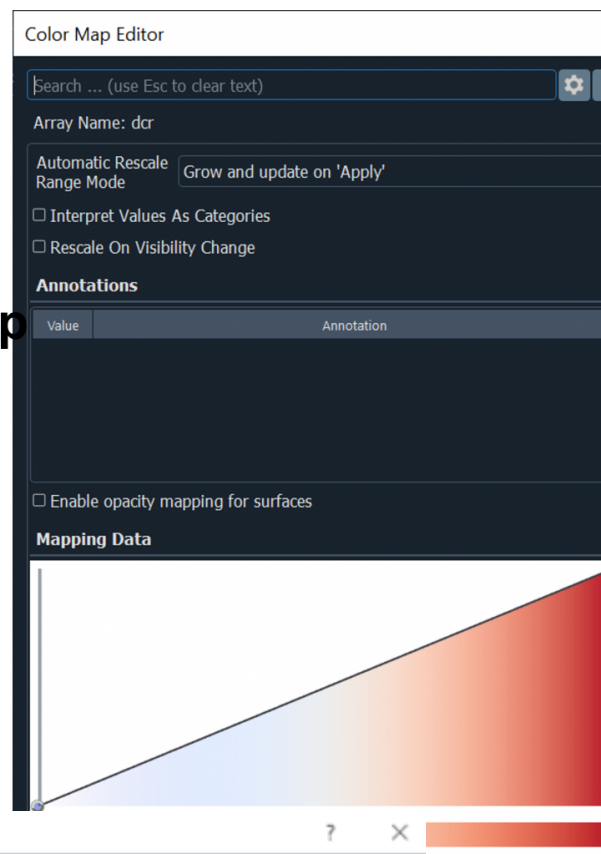


8. Choose Preset

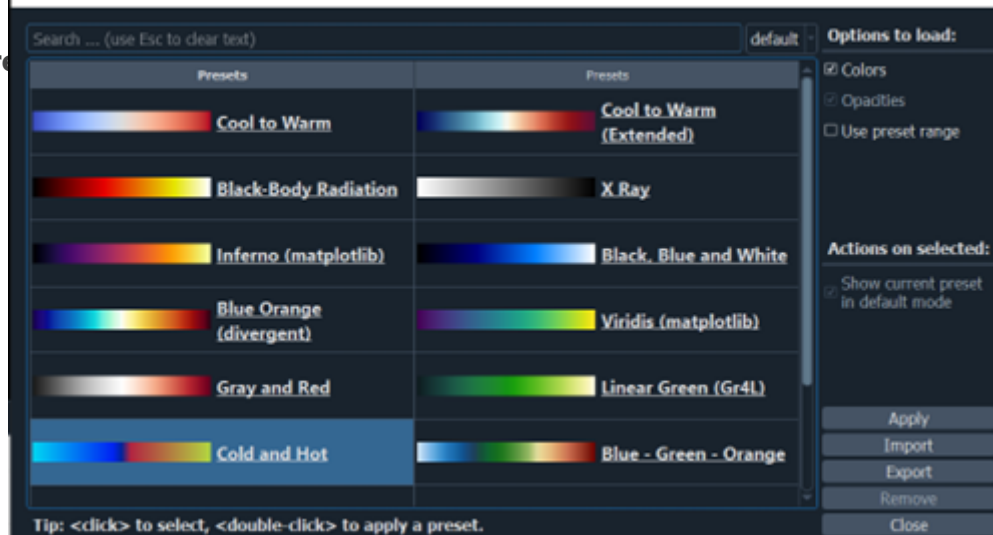


and click

apply

9. drag the color map **(A)** on the left side up **(B)**→ we want to discard dcr value from 0.32 to 0.36 we will set this to black color: set the left circle the lowest dcr value 0.32 (enter the value on Data and click enter) and the right circle to 0.36 → make another two circles 0.32 and 0.36 on the middle circle double click on the center of circle then select color: choose black, repeat for the left circle for example select cyan HTML#00e9e9, and for the right circle select orange HTML #ff8206

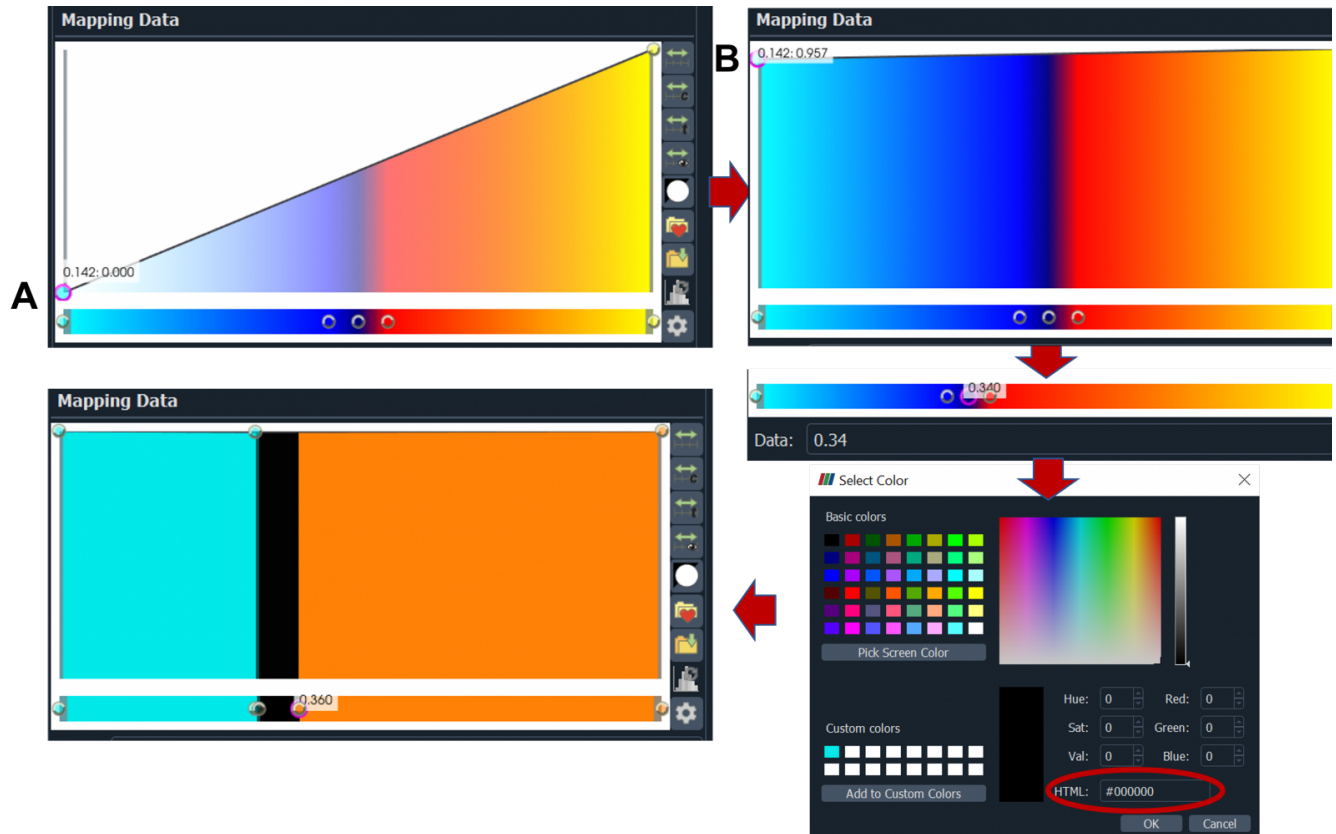

10. Signal with dcr value from 0.14 to 0.32 in cyan and dcr value of from 0.36 to 0.7 in orange

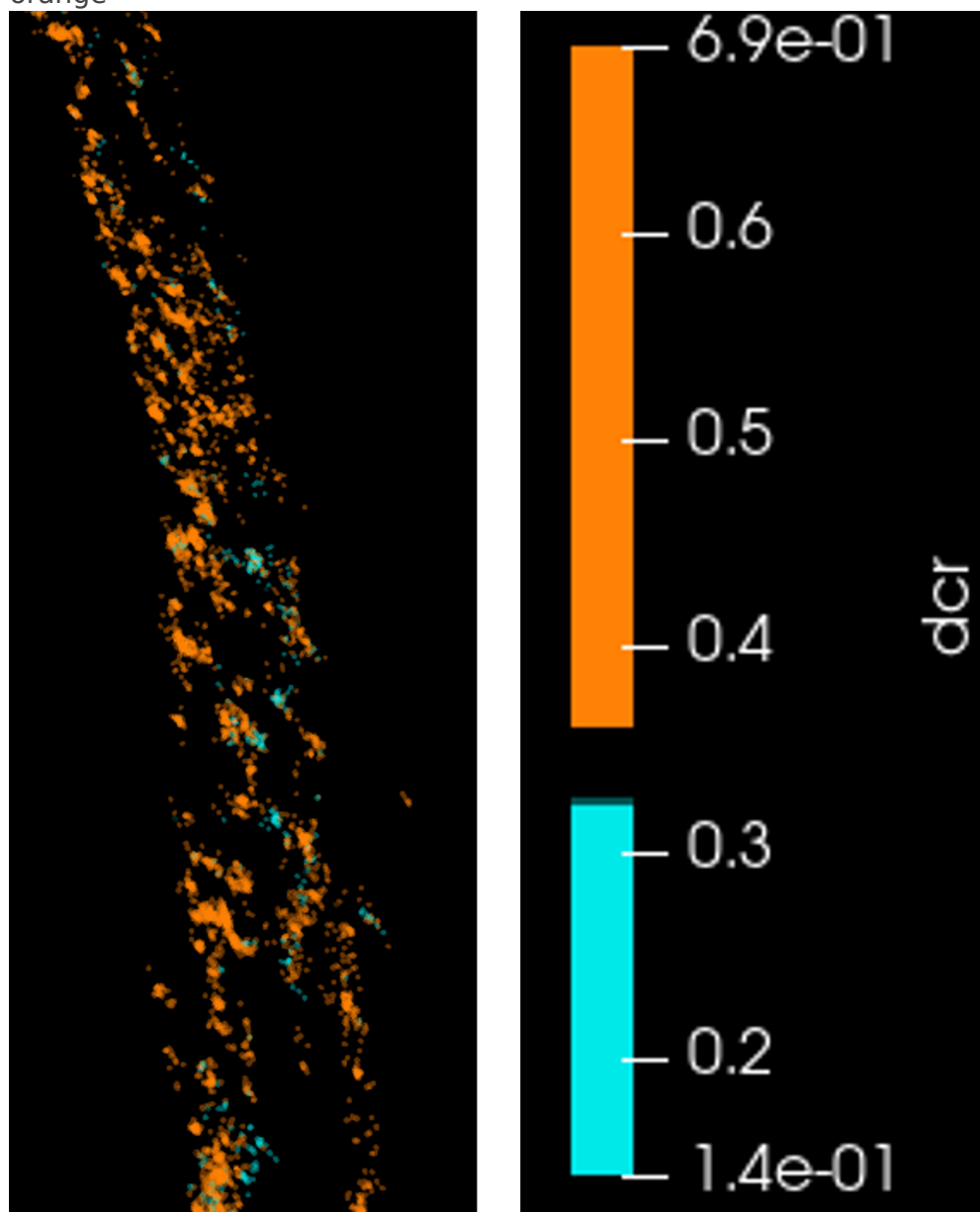

---

Revision #3

Created 28 July 2022 20:52:34 by Ni Putu Dewi Nurmalasari

Updated 17 January 2024 17:16:29 by Ni Putu Dewi Nurmalasari